



Multi-Array™ Immunodetection of Biowarfare Agents Using the Meso Scale Discovery™ Sector PR™

**Fred Lee¹, Deborah Menking¹, Joseph Hamel², Kishna Mangaya²,
Emily Myers¹, Michael Goode¹, and Alan Zulich¹**
Biosensors Team

¹U.S. Army, RDECOM, Edgewood Chemical Biological Center, APG, MD

²Geo-Centers, Inc., Aberdeen Proving Ground, MD



Report Documentation Page				Form Approved OMB No. 0704-0188	
Public reporting burden for the collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington VA 22202-4302. Respondents should be aware that notwithstanding any other provision of law, no person shall be subject to a penalty for failing to comply with a collection of information if it does not display a currently valid OMB control number.					
1. REPORT DATE 01 OCT 2005		2. REPORT TYPE N/A		3. DATES COVERED -	
4. TITLE AND SUBTITLE Multi-Array Immunodetection of BiowarfareAgents Using the Meso Scale Discovery Sector PR				5a. CONTRACT NUMBER	
				5b. GRANT NUMBER	
				5c. PROGRAM ELEMENT NUMBER	
6. AUTHOR(S)				5d. PROJECT NUMBER	
				5e. TASK NUMBER	
				5f. WORK UNIT NUMBER	
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) U.S. Army, RDECOM, Edgewood Chemical Biological Center, APG, MD 2Geo-Centers, Inc., Aberdeen Proving Ground, MD				8. PERFORMING ORGANIZATION REPORT NUMBER	
9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES)				10. SPONSOR/MONITOR'S ACRONYM(S)	
				11. SPONSOR/MONITOR'S REPORT NUMBER(S)	
12. DISTRIBUTION/AVAILABILITY STATEMENT Approved for public release, distribution unlimited					
13. SUPPLEMENTARY NOTES See also ADM001851, Proceedings of the 2003 Joint Service Scientific Conference on Chemical & Biological Defense Research, 17-20 November 2003. , The original document contains color images.					
14. ABSTRACT					
15. SUBJECT TERMS					
16. SECURITY CLASSIFICATION OF:			17. LIMITATION OF ABSTRACT UU	18. NUMBER OF PAGES 28	19a. NAME OF RESPONSIBLE PERSON
a. REPORT unclassified	b. ABSTRACT unclassified	c. THIS PAGE unclassified			



ECBC

ECBC Bioanalysis

- Current bioanalysis roadmap supports an integrated approach for detection and identification using multiple technologies
- We have combined state-of-the-art immunoassay analysis (ORIGEN[®] ECL) and nucleic acid analysis (SmartCycler[®] Taqman[™])
- Our current bioanalysis approach employs a tiered structure
 - Tier 1 : Presumptive ID by immunoassay
 - Tier 2 : Pathogen confirmation by PCR/Cell culture
Toxin confirmation by MS



ECBC

Immunoassay Challenges

- Front-end immunoanalysis requires development of sensitive & specific assays for presumptive screening of bioagents
 - Assay sensitivity cannot compromise assay specificity
 - Set assay cutoffs and LODs
 - Validate each assay for the matrix being analyzed
 - Cutoff must be above any sample matrix effects
 - Target sensitivities within the LD₅₀ and infective dosage ranges
 - Assay specificity has been an issue for some immunoassays
 - Need affinity-based antibodies directed to specific targets
 - Need to overcome sample matrix issues
 - Need to investigate alternative front-end sample clean up methodologies



ECBC

Immunoassay Challenges, cont.

- Identify technologies that feature simultaneous analyses of samples for multiple targets
 - Demonstrate mobility and simplicity for field use, urban monitoring, and mobile labs
 - Identify reliable, simple, rapid, and cost effective immunoassays
- Identify technologies that yield higher throughput analysis



Electrochemiluminescence

- ECL technology is elegant and yields exquisite sensitivity
- ECL measurement is simple and rapid
- Ruthenium tris (bipyridyl) $\text{Ru}(\text{bpy})_3^{2+}$ is a small molecule (< 1000 Da)
- Coupling chemistry is easy and straightforward
- $\text{Ru}(\text{bpy})_3^{2+}$ is an extremely stable compound



ECBC ECL FASTube History



IGEN ORIGIN



FY95-96

FY97-98

FY99

FY00

FY01-03

- R&D Evaluation COTS
- Optimize wet chemistry
- Publish results
- Develop one-step ECL assay
- Dry down reagents
- JFT assessment

FASTube validation

- U.S. Patent Application filed for method and device
- FASTube transition to TEU

FASTube technology transition to USAMRIID

- 2001 FASTube technology transitions to JPEO-CBD
- 2002 full-scale FASTube production by IGEN

- 2003 FASTube assays support Homeland Defense, TAML, TEU, OGAs, Iraqi Freedom

- 2003 FASTube assays offered in JPEO-CBD Reagent Catalog



Sector PR™ Pattern Array ECL Technology


- The Sector PR™ is a compact ECL reader designed by Meso Scale Discovery™ (MSD)
- The Sector PR™ is shoebox-sized and reads ECL signals from a 96 well plate in less than two minutes





ECBC

Multi-Array™ Technology: Printed Pattern Electrodes

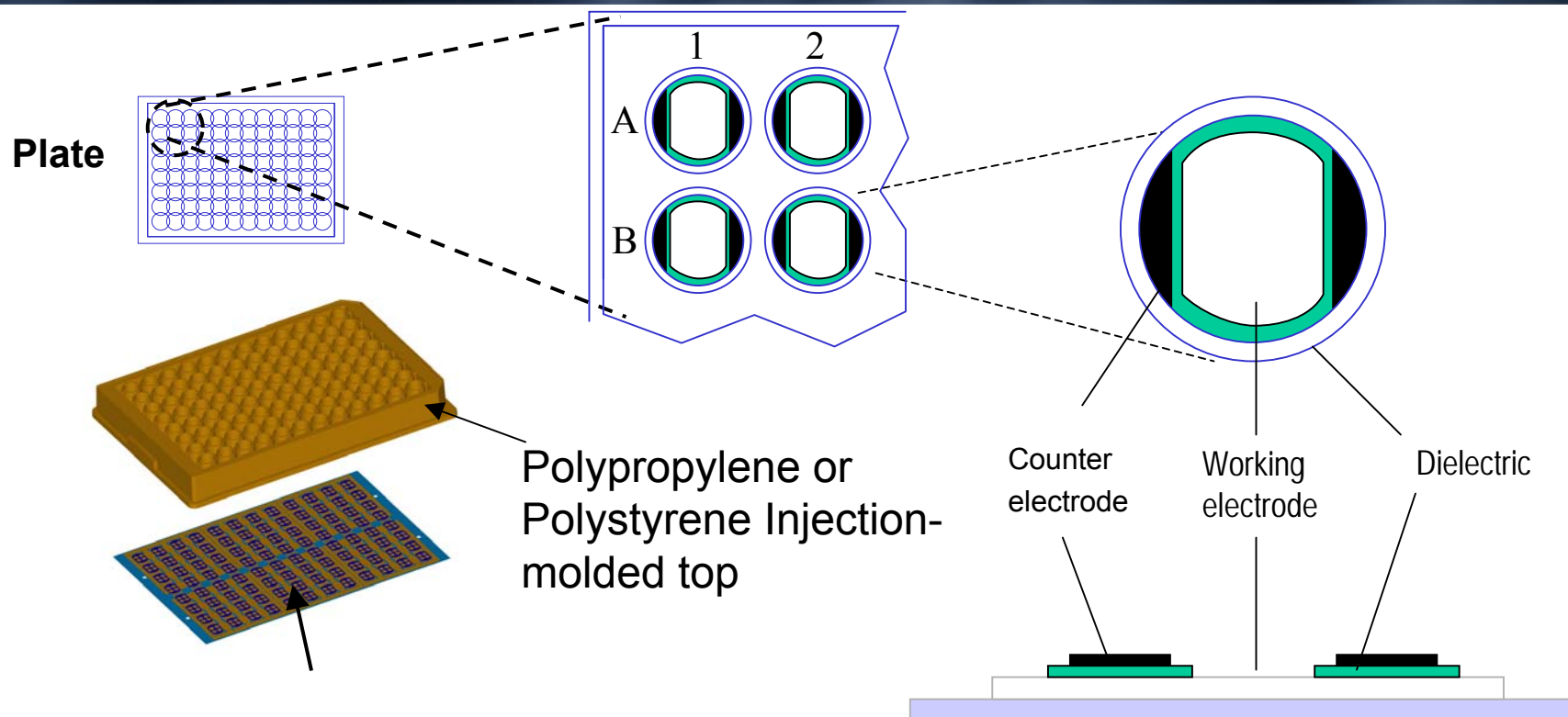

Screen Printed High Density Electrode Arrays

- Arrays are made from several patterned layers including layers of conducting carbon inks for the electrodes



ECBC

Objectives

- Evaluate Sector PRTM ECL assays:
 - Assay endpoint sensitivities and dynamic range
 - Intra- and Inter-assay variation
 - Assay specificity
 - Analysis of environmental samples (ENV)
 - Artificially spiked ENV samples
- Demonstrate proof-of-principle Multi-ArrayTM PR immunoassays
 - Staphylococcal enterotoxin B toxin
 - *Bacillus anthracis* spores
 - *Yersinia pestis*, F1
 - Venezuelan Equine Encephalitis
- Compare Sector PRTM results to the gold standard assays conducted on the ORIGEN[®] 1.5 ECL sensor



ECBC


Reaction chemistries

- Immobilize biotin conjugated capture antibody onto MSD streptavidin-coated 96 well carbon electrode plates
- Add ruthenium conjugated detector antibody and antigen to form an immunocomplex on the surface of the carbon electrode
- Apply voltage to the electrodes in the presence of the co-reactant tripropylamine (TPA)
- Eight silicon diodes sense the release of photons generated from the applied potential to the electrode, which is converted to the ECL signal



ECBC

ECL Detection

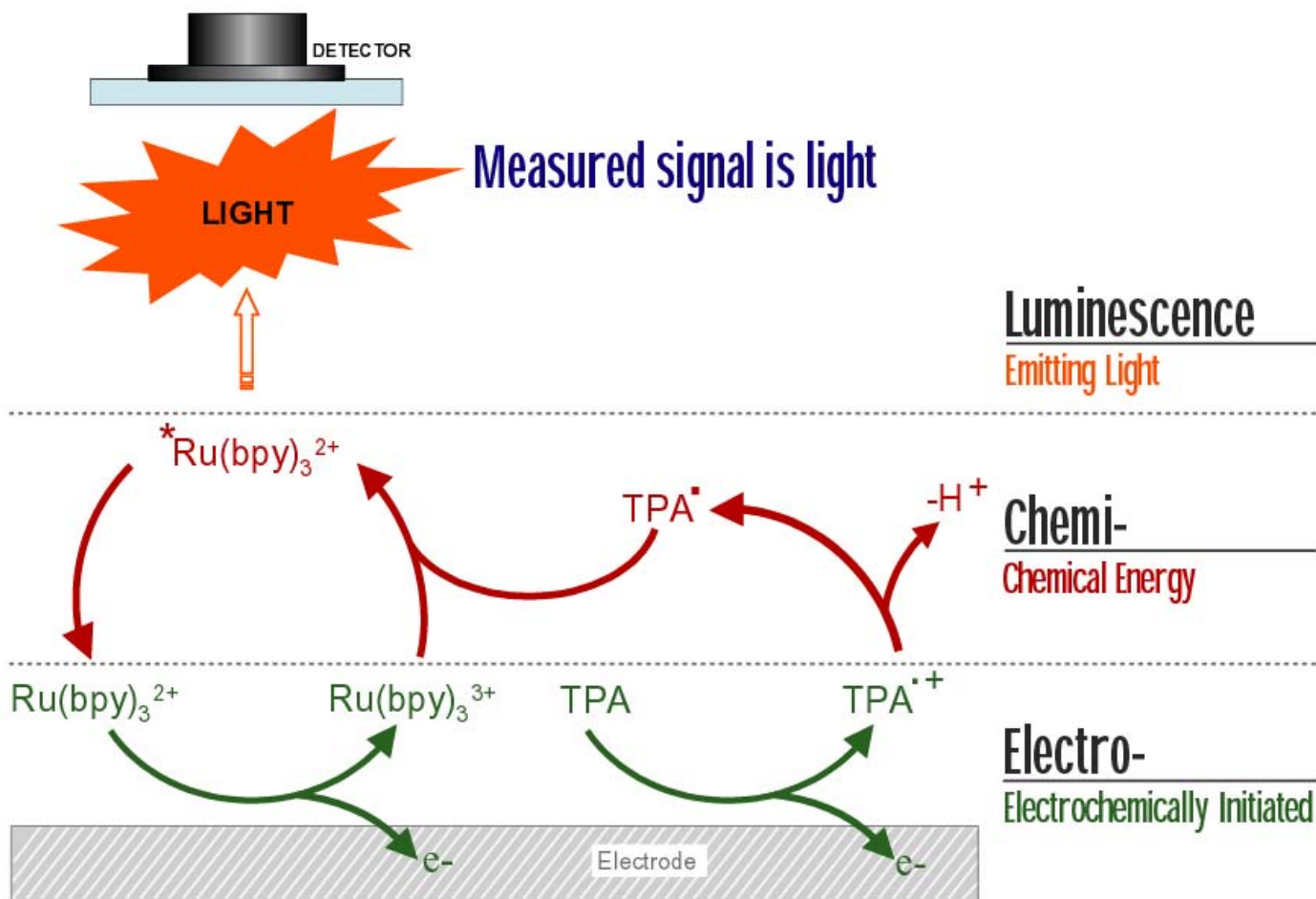

Meso Scale Discovery

A division of Meso Scale Biosciences, LLC



ECBC

Materials and Methods

- Antibodies:
 - Government furnished antibodies
 - Antibody conjugates
 - biotinylated capture antibodies
 - detector antibody labeled with ruthenium tris bipyridyl
- Format:
 - Mab anti-SEB capture / Pab rabbit anti-SEB detector
 - Mab anti-*B. anthracis* capture / Pab rabbit anti-*B. anthracis* detector
 - Mab anti-VEE capture / Mab anti-VEE detector
 - Mab anti-*Y. pestis* capture / Pab rabbit anti-*Y. pestis* detector



ECBC

Materials and Methods, cont.

Antigens:

- SEB toxin, Sigma, (0 -100 ng/ml)
- Gamma-inactivated *Bacillus anthracis*, Ames strain, AFIP, (0 - 1×10^6 cfu/ml)
- Gamma-inactivated VEE virus, TC-83 strain, USAMRIID, (0 - 1×10^7 pfu/ml)
- Gamma-inactivated *Yersinia pestis*, India strain, DPG, (0 - 1×10^6 cfu/ml)



ECBC

Materials and Methods, cont.

Assay Optimization:

- Capture antibodies matrixed from 30 -120 ng/test
- Detector antibodies were matrixed in the same manner as the capture
- The optimal capture/detector antibody concentration was determined through examination of endpoint sensitivity, CV, and dynamic range within the 15 -minute assay timeframe



Materials and Methods, cont.

Assay protocol:

- Coat MSD Streptavidin-Coated Multi-Array™ PR plate wells with 25 μ l of capture antibody
- Incubate at 4°C ON and wash three times in 0.01M PBS prior to analysis
- Add 25 μ l of detector antibody
- Add 100 μ l of the analyte to well and incubate for 15 min
- Wash plate three times with 0.01M PBS
- Add 150 μ l of 1X MSD Read Buffer to each of the wells
- Read plate on the Sector PR

Sector PR™ Washed Assay Data

Background + 3x STD:			927 ECL Units	
Lower Limit of Detection:			0.01 ng/ml	
SEB MSD Washed Assay				
60 ng biotin / 50 ng ruthenium per test				
ECL				
	Average	ECL		
	Raw	Specific		
ng/mL	Signal	Signal	% CV	S/N
0	859	-68	2.6	1.0
0.001	831	-96	3.5	1.0
0.005	891	-36	3.4	1.0
0.010	1043	116	2.5	1.2
0.10	3005	2078	2.3	3.5
1.0	23356	22429	1.2	27.2
10.0	222535	221608	3.7	259.1
100.0	554304	553377	15.8	645.3

Background + 3x STD:			704 ECL Units	
Lower Limit of Detection:			5.0E+04 cfu/ml	
<i>B. anthracis</i> MSD Washed Assay				
90 ng biotin / 100 ng ruthenium per test				
ECL				
	Average	ECL		
	Raw	Specific		
cfu/mL	Signal	Signal	% CV	S/N
0	593	-111	6.2	1.0
1.0E+01	563	-141	2.8	0.9
1.0E+02	599	-105	9.9	1.0
1.0E+03	586	-118	6.3	1.0
1.0E+04	651	-53	5.1	1.1
5.0E+04	741	37	2.0	1.3
1.0E+05	880	176	2.9	1.5
1.0E+06	2306	1602	4.3	3.9

Background + 3x STD:			802 ECL Units	
Lower Limit of Detection:			5.0E+06 pfu/ml	
VEE MSD Washed Assay				
30 ng biotin / 100 ng ruthenium per test				
ECL				
	Average	ECL		
	Raw	Specific		
pfu/mL	Signal	Signal	% CV	S/N
0	716	-86	4.0	1.0
1.0E+03	702	-100	4.8	1.0
1.0E+04	730	-72	10.5	1.0
1.0E+05	717	-85	4.4	1.0
5.0E+05	752	-49	5.3	1.1
1.0E+06	714	-88	1.7	1.0
5.0E+06	975	173	11.5	1.4
1.0E+07	1055	253	2.4	1.5

Background + 3x STD:			843 ECL Units	
Lower Limit of Detection:			1.0E+03 cfu/ml	
Y. pestis MSD Washed Assay				
60 ng biotin / 100 ng ruthenium per test				
ECL				
	Average	ECL		
	Raw	Specific		
cfu/mL	Signal	Signal	% CV	S/N
0	777	-66	2.8	1.0
1.0E+01	792	-51	3.1	1.0
1.0E+02	809	-34	4.2	1.0
1.0E+03	938	95	2.5	1.2
5.0E+03	1570	727	0.4	2.0
1.0E+04	2392	1549	1.5	3.1
1.0E+05	14349	13506	1.5	18.5
1.0E+06	58640	57797	3.5	75.4



Comparison of Platform Tradeoffs

Sensitivity		ORIGEN 1.5 [®]	Sector PR [™]
	<i>B. anthracis</i>	5 x 10 ³ cfu/ml	5 x 10 ⁴ cfu/ml
	<i>Y. pestis</i>	1 x 10 ² cfu/ml	1 x 10 ³ cfu/ml
	VEE	5 x 10 ⁵ pfu/ml	1 x 10 ⁶ pfu/ml
	SEB	5 pg/ml	10 pg/ml
Tradeoffs	Incubation Time	15 min	15 min
	Read Time	75 min	2 min
	# of Determinations	50	96
	Footprint	12'' x 12'' x 16''	8.6'' x 8.8'' x 16''
	Weight	47 lbs.	18 lbs.
	In-house cost per assay	\$5 – \$7	< \$2

Multi-Array™ PR Crossreactivity Results

Anthrax challenge

	SEB			<i>B. anthracis</i>			VEE			<i>Y. pestis</i>		
Bkgd	618	627	599	Flyer	640	689	625	676	806	935	917	997
5e04 cfu/ml	585	653	637	873	823	778	681	624	697	882	853	869
1e05 cfu/ml	632	587	557	921	964	924	667	606	674	875	888	862
1e06 cfu/ml	643	593	711	2436	2195	2256	658	622	642	892	894	896
Bkgd	612	625	623	598	612	682	592	601	619	767	738	777
0.010 ng/ml	775	796	789	651	699	691	700	670	687	893	847	913
1 ng/ml	22598	19526	20056	692	714	719	627	726	732	923	952	917
100 ng/ml	677958	658647	634486	699	783	804	767	732	754	890	933	958

SEB challenge



ECBC

Crossreactivity Results, cont.

VEE challenge

	SEB			<i>B. anthracis</i>			VEE			<i>Y. pestis</i>		
Bkgd	531	594	531	681	614	625	608	600	620	763	763	779
1e06 pfu/ml	585	593	598	673	675	671	669	667	685	770	785	770
5e06 pfu/ml	644	636	610	678	675	611	872	880	880	801	814	822
1e07 pfu/ml	607	572	523	679	642	631	984	1000	979	760	730	836
Bkgd	545	586	539	617	637	676	621	673	617	738	766	792
1e03 cfu/ml	557	605	550	664	656	646	605	618	651	1046	1089	1105
1e04 cfu/ml	567	608	597	697	669	698	679	676	670	3502	3552	3640
1e06 cfu/ml	508	573	575	645	648	629	658	659	610	111711	109926	117330

Y. pestis challenge



Environmental Sample Analysis

Multi-Array™ PR Plate Set Up

ENV sample	SEB			<i>B. anthracis</i>			VEE			<i>Y. pestis</i>		
Neg. Cntrl.	386	396	389	393	462	428	376	385	413	548	524	465
ENV 1	452	398	402	498	445	501	439	457	420	609	640	597
ENV 2	396	383	347	436	425	437	329	392	340	571	560	499
ENV 3	482	450	399	488	503	487	400	433	416	628	644	644
ENV 4	491	497	418	469	516	475	437	457	473	654	629	620
ENV 1 Spike	25316	25705	24561	1661	1673	1709	997	1059	1019	3065	2940	2897
ENV 3 Spike	23819	23209	22361	1579	1569	1482	827	901	826	2779	2728	2697
Pos. Cntrl.	23224	22839	22499	1642	1553	1466	1248	1297	1254	2490	2381	2484



ECBC

Environmental Sample Analysis Results

SEB

Day 1

Sample	ECL Raw Signal	S/N	Result
buffer negative	371	1.0	-
ENV-1	418	1.1	-
ENV-2	379	1.0	-
ENV-3	449	1.2	-
ENV-4	455	1.2	-
ENV-1 spiked	30181	81.3	+
ENV-3 spiked	28281	76.2	+
positive control	27099	73.0	+

Day 2

Sample	ECL Raw Signal	S/N	Result
buffer negative	390	1.0	-
ENV-1	417	1.1	-
ENV-2	375	1.0	-
ENV-3	444	1.1	-
ENV-4	469	1.2	-
ENV-1 spiked	25194	64.5	+
ENV-3 spiked	23130	59.3	+
positive control	22854	58.5	+

Day 3

Sample	ECL Raw Signal	S/N	Result
buffer negative	402	1.0	-
ENV-1	483	1.2	-
ENV-2	403	1.0	-
ENV-3	494	1.2	-
ENV-4	499	1.2	-
ENV-1 spiked	16942	42.1	+
ENV-3 spiked	16044	39.9	+
positive control	19064	47.4	+

B. anthracis

Day 1

Sample	ECL Raw Signal	S/N	Result
buffer negative	446	1.0	-
ENV-1	471	1.1	-
ENV-2	449	1.0	-
ENV-3	505	1.1	-
ENV-4	494	1.1	-
ENV-1 spiked	1623	3.6	+
ENV-3 spiked	1701	3.8	+
positive control	1758	3.9	+

Day 2

Sample	ECL Raw Signal	S/N	Result
buffer negative	428	1.0	-
ENV-1	481	1.1	-
ENV-2	433	1.0	-
ENV-3	493	1.2	-
ENV-4	487	1.1	-
ENV-1 spiked	1681	3.9	+
ENV-3 spiked	1543	3.6	+
positive control	1554	3.6	+

Day 3

Sample	ECL Raw Signal	S/N	Result
buffer negative	476	1.0	-
ENV-1	525	1.1	-
ENV-2	435	0.9	-
ENV-3	502	1.1	-
ENV-4	509	1.1	-
ENV-1 spiked	1621	3.4	+
ENV-3 spiked	1549	3.3	+
positive control	1529	3.2	+



Environmental Sample Analysis Results

VEE

Day 1

Sample	ECL Raw Signal	S/N	Result
buffer negative	391	1.0	-
ENV-1	420	1.1	-
ENV-2	403	1.0	-
ENV-3	437	1.1	-
ENV-4	420	1.1	-
ENV-1 spiked	929	2.4	+
ENV-3 spiked	836	2.1	+
positive control	1165	3.0	+

Day 2

Sample	ECL Raw Signal	S/N	Result
buffer negative	391	1.0	-
ENV-1	439	1.1	-
ENV-2	354	0.9	-
ENV-3	416	1.1	-
ENV-4	456	1.2	-
ENV-1 spiked	1025	2.6	+
ENV-3 spiked	851	2.2	+
positive control	1266	3.2	+

Day 3

Sample	ECL Raw Signal	S/N	Result
buffer negative	427	1.0	-
ENV-1	489	1.1	-
ENV-2	348	0.8	-
ENV-3	465	1.1	-
ENV-4	478	1.1	-
ENV-1 spiked	984	2.3	+
ENV-3 spiked	837	2.0	+
positive control	1457	3.4	+

Y. pestis

Day 1

Sample	ECL Raw Signal	S/N	Result
buffer negative	562	1.0	-
ENV-1	611	1.1	-
ENV-2	554	1.0	-
ENV-3	627	1.1	-
ENV-4	641	1.1	-
ENV-1 spiked	3040	5.4	+
ENV-3 spiked	2920	5.2	+
positive control	2693	4.8	+

Day 2

Sample	ECL Raw Signal	S/N	Result
buffer negative	512	1.0	-
ENV-1	615	1.2	-
ENV-2	543	1.1	-
ENV-3	639	1.2	-
ENV-4	634	1.2	-
ENV-1 spiked	2967	5.8	+
ENV-3 spiked	2735	5.3	+
positive control	2452	4.8	+

Day 3

Sample	ECL Raw Signal	S/N	Result
buffer negative	583	1.0	-
ENV-1	684	1.2	-
ENV-2	549	0.9	-
ENV-3	814	1.4	-
ENV-4	717	1.2	-
ENV-1 spiked	3859	6.6	+
ENV-3 spiked	3741	6.4	+
positive control	2935	5.0	+



Wrap Up

- We have developed four 15-minute washed assays on the Sector PR, a novel bioassay platform from Meso Scale Discovery
- Assays demonstrated excellent sensitivity and specificity
- Overall sensor performance was comparable to the ORIGEN[®] 1.5
- Exquisite CV's were achieved - below 5% in most cases
- Assay protocols have been developed for both washed and nonwashed formats. The nonwashed format is under development



ECBC

Wrap Up, cont.



- ENV samples were challenged with four antigens in Multi-Array fashion
- Cross-reactivity studies performed on the Multi-Array plates showed no evidence of nonspecific interaction
- The Sector PR assays are very sensitive, robust, and cost effective
- The Sector PR is shoe-box sized, is light and portable, and has no fluidics



ECBC

Future Objectives

Plate

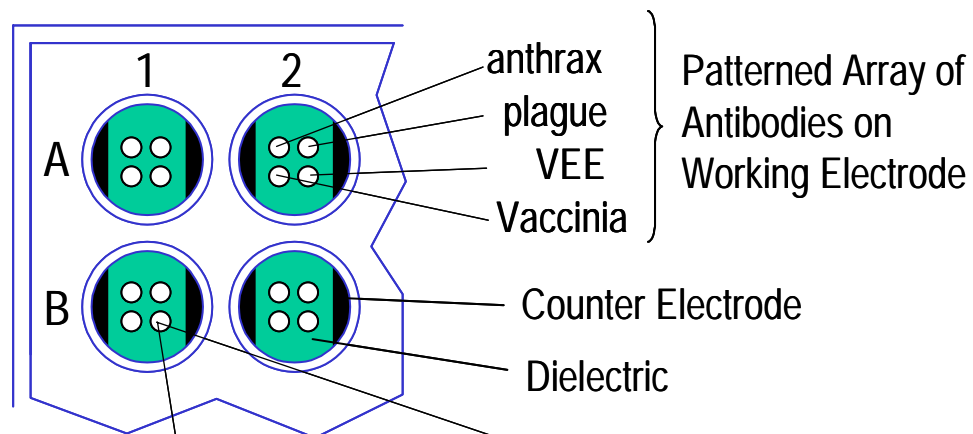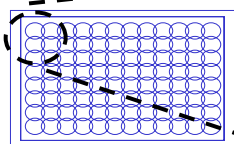


- **Multi-spot Pathogen Test**

- B. anthracis*
- VEE
- Y. pestis*
- Vaccinia

- **Multi-spot Toxin Test**

- Ricin
- SEB
- C. botulinum*

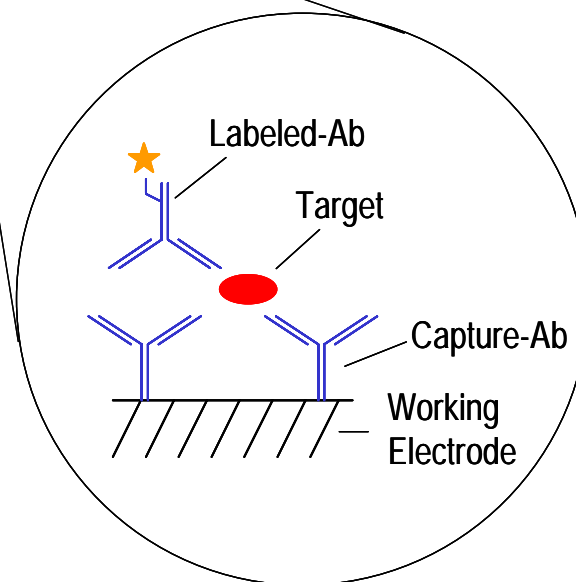



Future Objectives, cont.

- Investigate feasibility of dried-down reagents in MSD 96-well plates using StabilCoat® buffer and lyophilization buffer
- Proceed with assay development of additional biowarfare agents
- Refine wash and nonwashed assay formats
- Continue to evaluate MSD next generation bioanalysis platforms
- Evaluate sensor potential for portability, ease of use, sensitivity, and speed



ECBC

Technology Watch

Instrument



**Ruggedized
Sector PR**




Strip Reader




Cartridge Reader

Consumables

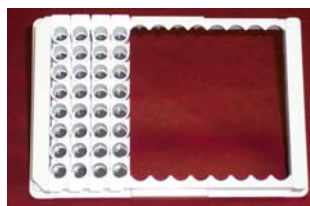

**96-Well Plate,
4-Spot Plate,
Multiple Strips
of Wells.
Dried reagents
in wells**

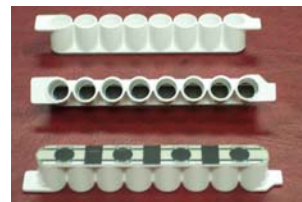

**Single Strip
of Wells.
Dried
reagents
in wells**



**Cartridge
contains
all consumables**



ECBC

Acknowledgements

Funding was provided through an ILIR grant from ECBC Research & Technology Directorate.

ECBC Biosensors Team



Meso-Scale Discovery™



Mr. Vittal Vasista

Dr. Eli Glezer

Dr. George Sigal



Meso Scale Discovery®
A Division of Meso Scale Biosciences, LLC.